TERMOSIGILLATRICE ELETTRONICA
ELECTRONIC THERMOSEALING MACHINE
ELEKTRONISCHES FOLIENSCHWEIßGERÄT
THERMOSOUDEUSE ELECTRONIQUE
TERMOSELLADORA ELECTRÓNICA
APARELHO ELECTRÓNICO DE SELAGEM TÉRMICA

EUROSEAL® 2001 PLUS

MANUALE DI ISTRUZIONE E MANUTENZIONE
INSTRUCTION MANUAL AND MAINTENANCE
GEBRAUCHS UND INSTANDHALTUNG HANDBUCH
MANUEL D'INSTRUCTIONS
MANUAL DE INSTRUCCIONES
MANUAL DE INSTRUÇÃO E MANUTENÇÃO

Italiano

English

Deutsch

Français

Español

Português

Copyright © 2001 Euronda S.p.A.

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Dear Doctor,

thank you for the faith that you have given us by acquiring our thermosealing machine **EUROSEAL® 2001 PLUS**.

We wish to inform you that our firm is at your complete disposition to furnish information and to answer any questions concerning the operation and use of this equipment.

We remind you that for the correct usage of this machine it is necessary to first read this manual. The **EUROSEAL® 2001 PLUS** corresponds to all of the relative safety regulations and presents no danger to the operator when used according to the instructions.

While wishing you success in your work, we remind you that the reproduction of this manual is forbidden and that the technical characteristics could be changed without warning due to our continued technical research.

EURONDA S.p.A.

• GENERAL INDEX •

• INTRODUCTION

GENERAL NOTES ON DELIVERY	P. 3
GENERAL INSTRUCTIONS	P. 3
EXPLANATION OF SYMBOLS	P. 4
WARRANTY INFORMATION	P. 4

• PART 1

DIMENSIONS OF THE MACHINE AND OF THE PACKAGE	P. 5
--	------

• PART 2

INSTALLATION	P. 6
INSTRUCTION FOR WALL MOUNTING	P. 6
SAFETY	P. 7
ELECTRICAL CONNECTION	P. 8

• PART 3

DESCRIPTION OF THE THERMOSEALING MACHINE	P. 9
CHARACTERISTICS OF THE THERMOSEALER	P. 9
STERILISATION ROLLS AND BAGS	P. 10

• PART 4

INSTRUCTIONS FOR USE	P. 11
----------------------	-------

• PART 5

CLEANING INSTRUCTIONS	P. 13
SPECIAL MAINTENANCE AND WARNINGS	P. 13
HOW TO REPLACE THE RESISTANCE SUPPORT	P. 13
HOW TO REPLACE THE SILICON RUBBER BAND	P. 14
HOW TO REPLACE THE ELECTRONIC CIRCUIT BOARD	P. 14
HOW TO REPLACE THE CUTTING EDGE	P. 14
PROBLEM SOLVING GUIDE	P. 15

• PART 6

ELECTRICAL DIAGRAM	P. 16
DIAGRAM OF THE MACHINE COMPONENTS	P. 17
COMPONENT LIST	P. 18
SCRAPPING INSTRUCTIONS	P. 19

GENERAL NOTES ON DELIVERY

Upon receipt of the equipment, please control the condition of the package and keep it for possible shipments.

Open the package and check that:

- the equipment corresponds to the specific requirements;
- there is no evident damage.

In case of damage or missing parts, inform immediately either the shipping agent, the dental depot or Euronda S.p.A.

All of the dimensions listed in this manual are not binding.

The sketches and any other consigned documents regarding this machine are property of Euronda S.p.A., which reserves all rights. They can not be made available to any third parties.

It is also forbidden to reproduce in whole or in parts either the instruction text and the illustrations.

GENERAL INSTRUCTIONS

Included in this manual are instructions for:

- correct installation;
- the safe and efficient use of the apparatus;
- routine maintenance.

The unit must be used in compliance with the procedures contained in the manual and never for reasons other than those specified therein.

The **operator** is the person who physically uses the apparatus for the purposes for which it has been designed.

The **responsible authority** is the person or group responsible for its use, for the ordinary maintenance of the apparatus and for the training of the operator.

The responsible authority is to be held legally responsible for the performance of the unit as far as it is a direct result of its installation, operation and use.

The manufacturer cannot be held responsible for any breakage, damage or malfunctioning of the unit if the machine has not been used correctly or not adequately maintained.

Conditions of use: this manual is an integral part of the product and must be kept near the unit in order for it to be easily and quickly consulted. The unit is for professional use only and only qualified persons may use it. The unit is to be used only for the purpose for which it was designed.

**EUROSEAL® 2001 PLUS: is a machine for thermosealing and cutting of sterilisation rolls made of polypropylene paper.
It does not work with other materials.**

EXPLANATION OF SYMBOLS

Safety symbols present on the apparatus:

	WARNING: HIGH TEMPERATURE
---	----------------------------------

Some symbols are displayed on the label of the unit. Their meaning is illustrated, as follows:

SYMBOL	DESCRIPTION
SN	“SERIAL NUMBER” The symbol must be accompanied by the manufacturer’s serial number. The serial number must be adjacent to the symbol.
	“PRODUCTION DATE” The symbol must be accompanied by the year. The year must be composed of 4 digits
	“WARNING, READ THE INSTRUCTION MANUAL”

In the manual there are, in addition, symbols adjacent to a description or a note, etc. These symbols are intended to draw the attention of the reader to a particular note or explanation. Their meanings are explained below.

SYMBOL	DESCRIPTION
	IMPORTANT SAFETY INFORMATION This symbol is used to draw the attention of the reader to notes of particular importance for the safety of the user.

GUARANTEE INFORMATION


Euronda guarantees the quality of its equipment, if used in accordance with the instructions supplied in this manual, and according to the conditions printed on the GUARANTEE CERTIFICATE.

WARNING: the CUSTOMER USER must fill in all parts of the DETACHABLE SLIP of the guarantee certificate and send it to EURONDA S.p.A.


The guarantee is effective from the delivery date of the unit to the customer, which is confirmed by the return of the guarantee slip correctly filled in and signed. In the case of dispute, the date indicated on the purchase invoice, which shows the serial number of the unit, will be considered valid.

NOTE: Please retain the original packaging and use it each time the unit needs to be transported. The use of different packaging could damage the product during transportation.

PART 1 DIMENSIONS OF THE MACHINE AND OF THE PACKAGE

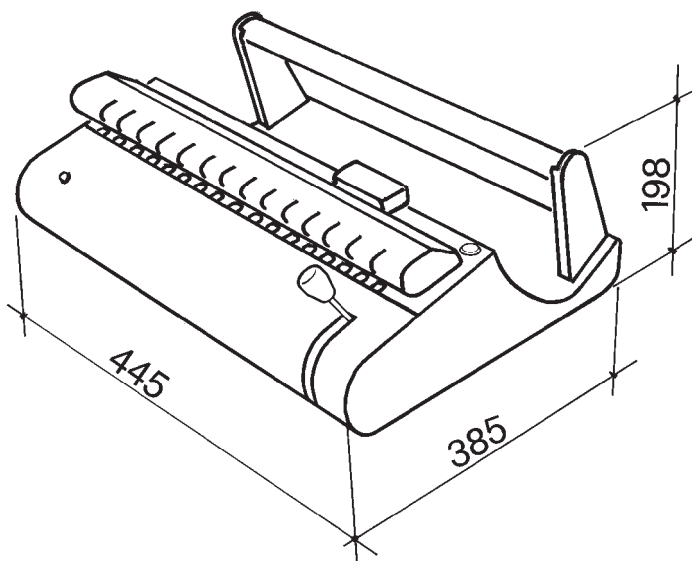Dimensions of the electronic thermosealing machine Euroseal® 2001 PLUS

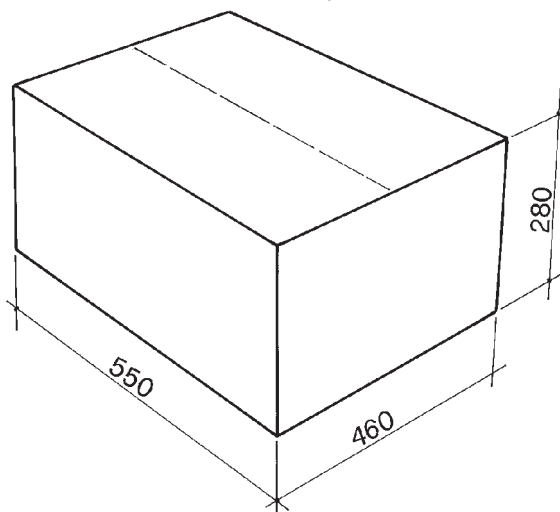Overall dimension of the package

PART 2

INSTALLATION


The thermosealing machine EUROSEAL® 2001 PLUS has been tried and tested in our factory.

Remove the packing and install the apparatus using the following instructions:

1. Install the machine on a plain flat surface with a remaining space of at least 50-80 mm. on all sides.
2. Do not install the machine close to a sink or similar to avoid contact with water and possible electrical short circuits.
3. Keep:± far from heat sources
4. Install the machine in a ventilated area.
5. Ensure that the electrical cable is not folded up and freely reaches the electrical socket.
6. Insert the plug and the machine is ready for use.

Upon request the machine is able to work also with 110V, 50/60 Hz.

INSTRUCTIONS FOR WALL MOUNTING

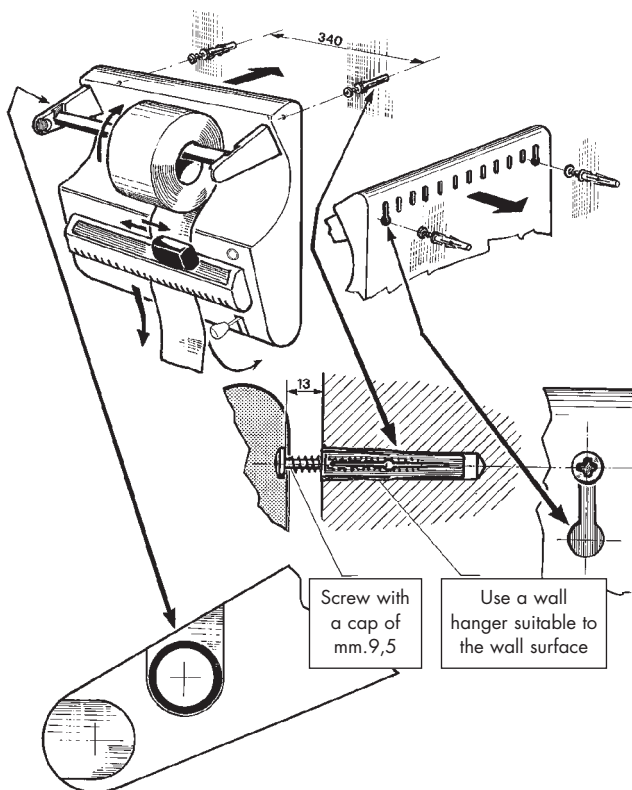SAFETY

Before using your electronic thermosealing machine EUROSEAL® 2001 PLUS, please read carefully this safety information. Your lack of observance could cause accidents or damage to the machine.

NEVER PLACE YOUR FINGERS NEAR THE SEALING AREA WHEN THE MACHINE IS OPERATING
NEVER PLACE YOUR HANDS UNDER THE SEALER WHEN YOU PRESS THE SEALING LEVER
NEVER PLACE YOUR HANDS NEAR THE SHARP EDGE OF THE CUTTING BLADE.

KEEP THE MACHINE IN A CLEAN AND DRY AREA.

NEVER DISMANTLE THE Euroseal® 2001 PLUS. THE HIGH INTERNAL VOLTAGE IS DANGEROUS

If it is not possible to disconnect the electricity or if it is too far or not visible to the person who is making the repairs, attach the sign "Work in Process" on the switch after having turned it "OFF".

CLEAN THE MACHINE WITH A DAMP CLOTH, AFTER ENSURING THAT THE CONNECTION CABLE HAS BEEN UNPLUGGED (before reusing the machine make sure that the machine is completely dry).

MAKE SURE THAT THE WIRING HAS A GROUND CABLE WHICH CONFORMS TO ALL OF THE SAFETY NORMS OF THE COUNTRY OF INSTALLATION.

NEVER REMOVE THE LABEL ON THE MACHINE

ENSURE THAT THE MACHINE IS CONNECTED TO THE CORRECT ELECTRICAL VOLTAGE

USE EXCLUSIVELY THE ORIGINAL REFILLS

If the machine doesn't WORK properly refer to the selection in the manual "Solution to Problems" for the possible cause. If you need subsequent information or the machine have to be repaired contact the specialised personnel of the supplier or the technical office of Euronda S.p.A.

WARNING

- The machine conforms to the safety regulations established by the Standards Institute and has a two prong plug which ensures that the apparatus is shock resistant.
- The device has a double isolation. Only the resistance needs to be enclosed in a base which is sufficient for conducts electricity..
- The most fundamental safety requirements are: a sturdy base and an electrical voltage capacity which is written on the tag. In case of doubt have it verified by qualified personnel.

**BY NOT OBSERVING THE REQUIREMENTS LISTED EURONDA S.p.A.
FORFEITS ALL RESPONSABILITY for its misuse.**

ELECTRICAL CONNECTION

Check that the electrical current indicated on the tag on the rear of the machine, corresponds to the voltage of the place of installation.

The machine must be connected to an electrical system with adequate voltage which corresponds to all standards laws of the country of installation.

ATTENTION: Don't fold the electrical cable and don't lean objects on it.

PART 3

DESCRIPTION

Euroseal® 2001 PLUS is the newest electronic thermostealing machine from Euronda, with a continuous cycle and constant temperature which has a guaranteed sealing size of 12 mm., and is available in both 220V or 110 voltage. Euroseal® 2001 PLUS is a machine with a modern, innovative and ergonomic design, easy to use and to clean, with rounded contours and fitted with a roll holder which has an overall dimension that measures the length of the machine. The sterilisation rolls of paper/polypropylene, available in 100 or 200 meters, are positioned inside the machine. The new applied technology minimises the energy consumption by creating shorter sealing times and eliminating the need of cooling times between every use. It also increases the reliability of the machine, excluding the periodical servicing and repairs. Euroseal® 2001 PLUS is built in accordance with the quality standards set by the European Community.

CHARACTERISTIC OF THE THERMOSEALER

Supply voltage	230 V	110 V
Number of phases	1	1
Supply frequency	50 Hz	50/60 Hz
Power input	100 W	100 W
Isolation class	1	1
Fuses	N° 2 - 5 x 20 - T 4 A	N° 2 - 5 x 20 - F 10 A
Range of environmental conditions in which the unit was designed to operate	<ul style="list-style-type: none"> - Indoor use - Altitude up to 2000 m - Temperature: +5 ÷ +40°C - Max relative humidity 85% - Max variation in mains voltage: ±10% - Installation category (overvoltage category) II - Pollution level 1 	
Switch	Single-pole	
Weight of the machine	7,4 Kg	
Sealing band	12 mm	
Type of sealing	constant temperature.	
Sealing area width	310 mm	
Power unit	unique and autonomous, available for use with either 220V or 110 voltage. Which is not subject to short circuiting. It is easy to repair (in the rare case of breakage), and is built according to EC safety standards.	
Roll hanger	has a hidden support	
Energy consumption	It is from 30% to 50% less than the standard impulse sealing mac	
Cooling time	not needed, between one sealing and another	
Positioning of the tape during the sealing	During the sealing there is no difference between placing the sealing tape with the film above or below	
Cutting edge	practical, safe, functional and easy to use	
Sealing lever	ergonomic, handle with a light and delicate movement	
Design	ergonomic and rounded	
Colour	RAL 9002, powder painted	
Routine maintenance	not necessary	

STERILISATION ROLLS AND BAGS

EUROSTERIL® sterilisation rolls and bags are obtained from high grain paper and polypropylene film. They ensure the hygienic storage of instruments both in envelopes and packages originated by a thermosealing machine.

The correct completion of the sterilisation process is signalled by a special indicator, when the colour of the strip on the side of the paper changes.

The roll consumption can be regulated based on the actual dimensions of the objects which are being sealed reducing, in this way, any waste.

The EUROSTERIL® rolls are available in the following sizes:

- Rolls of 100 meters long, in the following width: 5.5 - 7.5 - 10 cm.
- Rolls of 200 meters long, in the following width: 5.5 - 7.5 - 10 - 15 - 20 - 25 - 30 cm.
- Folding rolls of 100 meters, width 20 x 5 cms.

The EUROSTERIL® bags are available in the following sizes:

- 55x250 mm
- 75x250 mm
- 100x250 mm

Rolls and bags are produced according to the European safety standards related to the materials used in the process.


PART 4

INSTRUCTIONS FOR USE

1. Connect the machine the electricity source using the electrical cord.


Press the lateral switch POWER - **A** -, located on the top of the unit.

The light - **B** -, found in the front of the machine will turn green when the optimal sealing temperature has been reached.

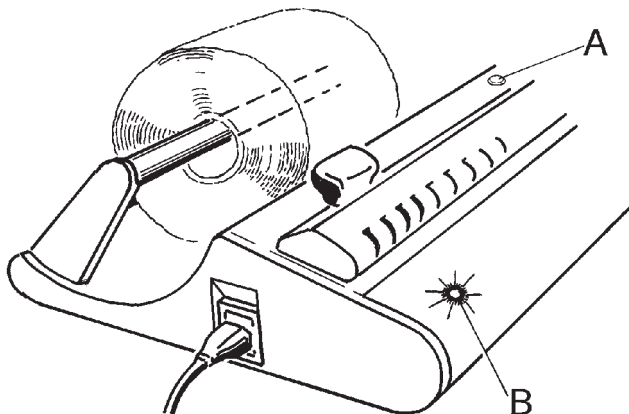
2. Fit the sterilisation roll under the holder -D-.

The roll can be inserted with the paper foil positioned either rightside up or upside down.

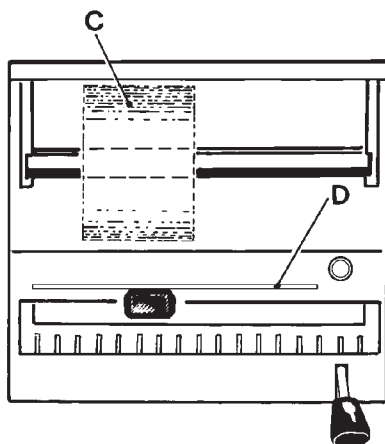3. Lower the sealing lever. The lever will pass from position -G- to position -H-.

When the sealing procedure is completed the light - I - will turn red and the signal bell will ring.

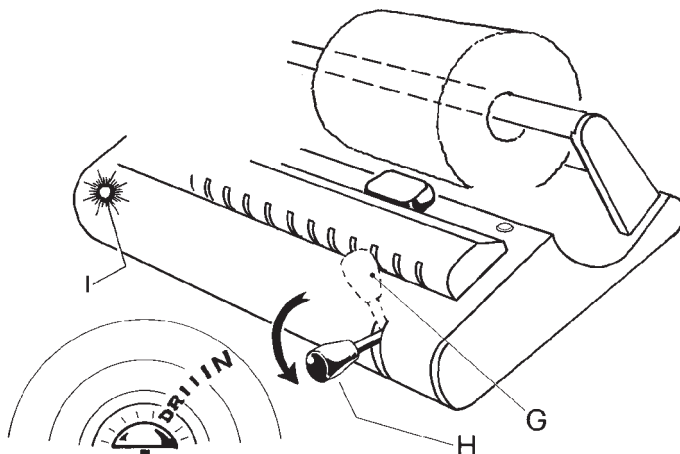
4. Raise the sealing lever passing from position -L- to position -M-.

When the lever is raised the light will turn green again and the signal bell will stop ringing.

The machine is ready for the next sealing without having to wait for a cooling period.

If the operator doesn't return the lever to its start position immediately after the sealing procedure is completed, the power unit will stop itself automatically (so the tape does not burn). The bell however will continue to ring. As soon as the lever is returned to its start position - G - the power unit will automatically switch on again. It will take several seconds to restore the optimal sealing temperature.

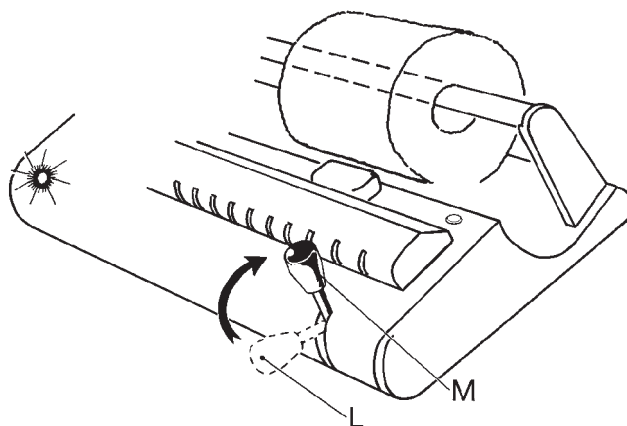PART 5

CLEANING INSTRUCTIONS

- Clean the exterior of the machine with a damp cloth, after making certain that the machine is completely dry. Before starting to work make sure that the machine is completely dry.
- During the cleaning of the machine make attention to the cutting edge position.

SPECIAL MAINTENANCE AND WARNINGS

- The machine does not require routine maintenance
- For every special maintenance operation refer to the specifications section in PART 2.
- **All operation must be effected by a specialised technician.**

HOW TO REPLACE THE POWER UNIT, SILICON RUBBER BAND AND OF THE ELECTRONIC CIRCUIT BOARD


RESISTANCE:

Remove the two plugs - **A** - from the knife protection casing and then loosen the two screws - **I** -; raise the whole unit - **G** - with the Belleville washers - **B** -.

Lower the lever - **C** -, as shown by the arrow, then remove the screw - **U** - from the resistance insulating block - **S** - and withdraw it by pulling it towards the left and lifting it until the connector - **E** - comes out from the resistance support.

It is now possible to disconnect the ground wire from machine connection - **D** - as well as the connector - **E** -.

To re-assemble, proceed inversely paying particular attention when inserting the wires into their seats. Re-fit the Belleville washers - **B** - before securing the unit - **G** - with the screws.

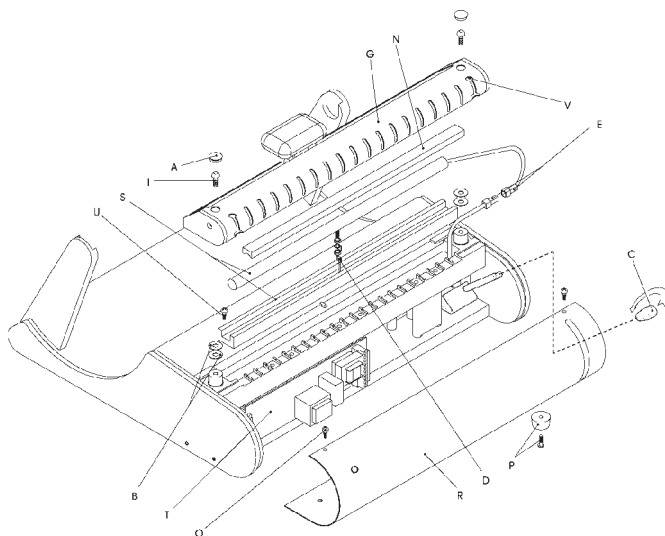SILICON RUBBER BAND:

Disassemble the unit - **G** - as described above and remove the silicone rubber - **N** -. Insert the new rubber by pressing it in.

ELECTRONIC CIRCUIT BOARD:


Lower the lever and unscrew the grip of pin - **C** -, remove the front carter - **R** - taking away the screws - **O** - acting with the screwdriver into the hole - **V** -. Remove also the two pins - **P** - under the machine.

Take away the screws - **T** - and the circuit card.

Once replaced the defective parts, reassemble the machine following the inverted order.

REPLACING OF THE CUTTING EDGE:

Remove the screw - **A** -, from one of the two sides - **B** - of the higher carrier - **C** - that simply slips off from the cuttinghandle - **E** -. Being careful not to cut yourself, slip off from the cutting bar - **D** -, the handle - **E** - and remove the two screws - **F** - and you can replace the blade - **G** -. Once the substitution is complete reassemble the machine in the reverse order.

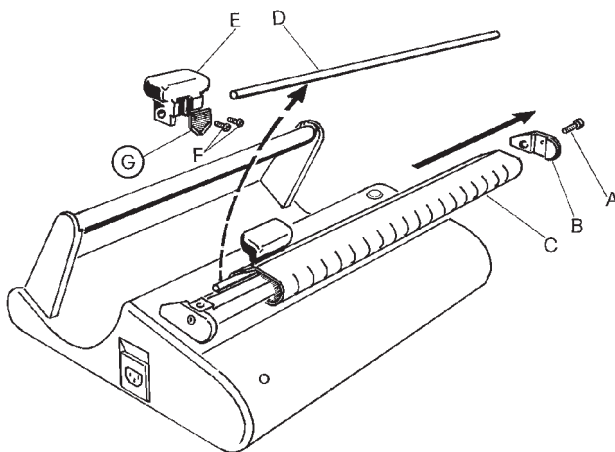PROBLEM SOLVING GUIDE


PROBLEM	CAUSES	SOLUTION
The switch light remains off.	No voltage	No electrical current
	Protection fuses on the side of the cable	Check for possible short circuit or breakage, if so replace with fuses of the same value
	Switch light broken	Replace the switch fuses set
The panel light doesn't turn on after the machine has been on for several minutes	Panel broken	Replace the electronic board
	The light doesn't function	Replace the light
	The power unit is not heating	Replace the power unit
The sealing process occurs normally, but the light remains green	Micro Switch broken	Replace the micro switch
The unit heats, but sealing process is not sufficient	Sealing pressure is not sufficient	Check the screws which support the silicon rubber piece
	Silicon rubber ruined	Replace the silicon rubber piece.

ATTENTION: BEFORE EVERY MAINTENANCE INTERVENTION TURN OFF THE VOLTAGE ALLOW A SPECIALISED TECHNICIAN TO PERFORM THE CHECK


DANGER: HIGH INTERNAL VOLTAGE

PART 6 ELECTRICAL DIAGRAM

English

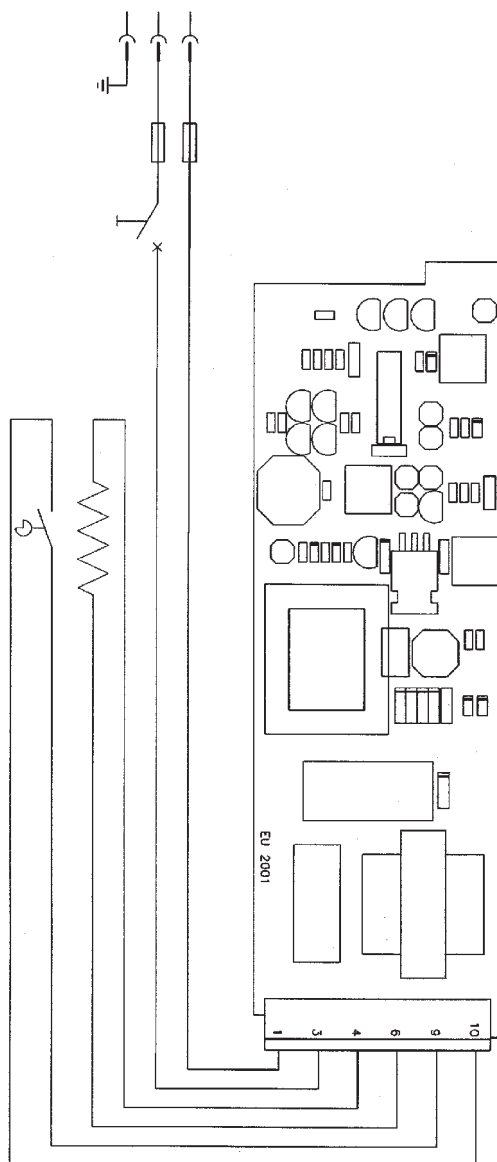COMPONENTS OF THE MACHINE

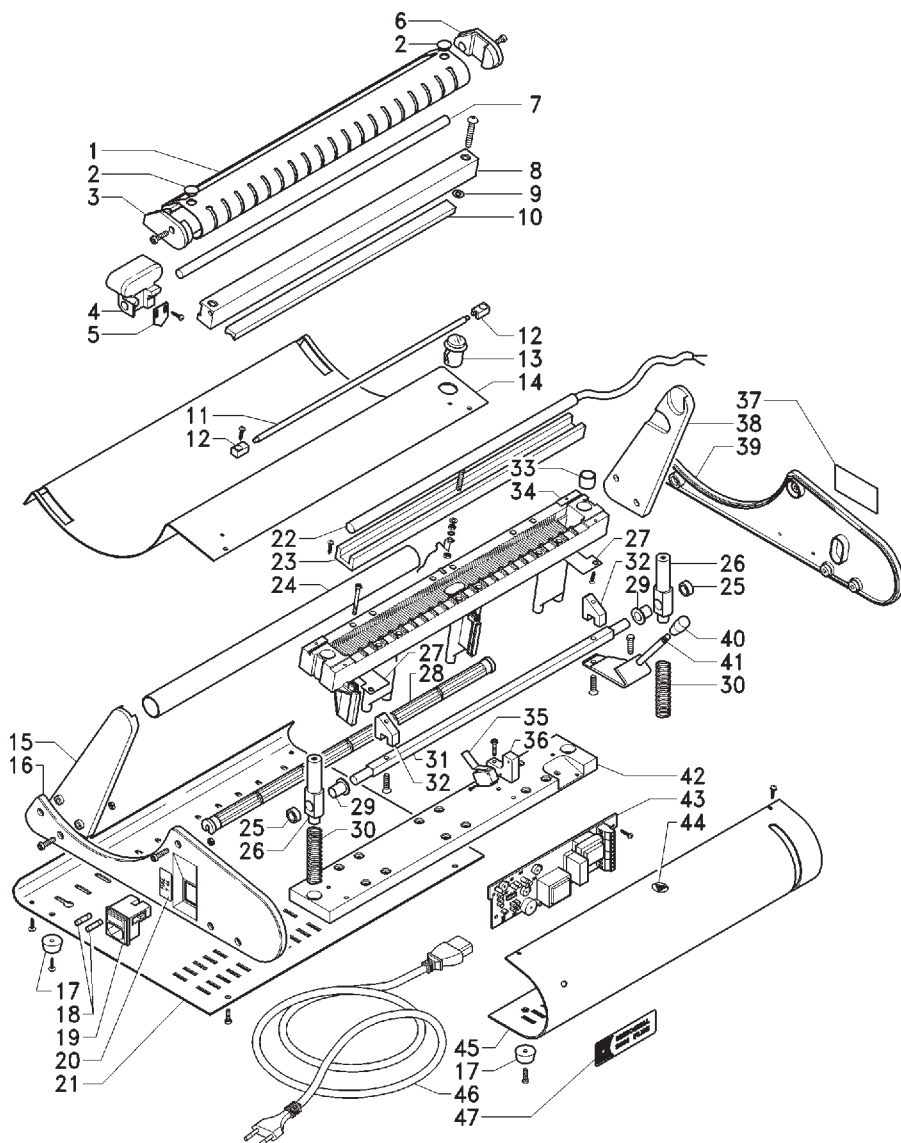COMPONENTS LIST

To order spare parts, please refer to the specific code.

NUMBER	QUANTITY	DESCRIPTION	CODE
1	1	Painted sheet metal knife guard	930125
2	2	Black rubber plug mm10	940214
3	1	LH side of upper guard. RAL 9002	931130
4	1	Black nylon knife guide	930128
5	1	Knife ric.X1052A 3-12 SCU	962201
6	1	RH side of upper guard. RAL 9002	931129
7	1	Stainless steel knife guide rod	930127
8	1	Anodised aluminium knife guide plate	930136
9	4	Belleville washer 12x5x0.6	962121
10	1	Silicon rubber contour 14x8 I.323	953212
11	1	Paper roller	930140
12	2	Roller support	930142
13	1	Switch	956085
14	1	Painted sheet metal upper guard	930185
15	1	Nylon LH roll holder support RAL9002	931113
16	1	Nylon LH side RAL9002	931002
17	4	Black rubber foot diam. 21	972138
18	2	Mod. A 230V: Fuse 5x20 T4A Mod. A 110V: Fuse 5x20 F10A	956117 956128
19	1	Socket with fuse holder	942039
20	1	Fuse label 38x18	948153
21	1	Painted sheet metal base	930116
22	1	Mod. A 230V: Heater 12x335 100W Mod. A 110V: Heater 12x335 100W	930201 930203
23	1	Glass fibre heater insulation 42	931117
24	1	Black PVC tube 25-22x390	930114
25	2	Nylon bush d.15x6	930131
26	2	Guide pin. 20x83	931120
27	2	Steel plate for cam 50x30	931139
28	1	Rod for sides RAL9002	930145
29	2	Self lubric. bush mb.15 20	936017
30	2	Spring diam. 15x58 wire 1.2	962153
31	1	Rod support for cams	930223
32	2	Cam in arnite	932123
33	2	Self lubric. bush fl. 10 15	936015
34	1	Heater support assembly	930147
35	1	Microswitch lever matsushita AH71525	851003
36	1	Nylon microswitch support RAL9002	930141
37	1	Print label 33x63 D.T. CE	948103
38	1	Nylon RH roll holder support RAL9002	931112
39	1	Nylon RH side RAL9002	931001
40	1	Knob I10825 TM601	995255
41	1	Steel lever 35x5 and pin in avp diam. 8x69	930138
42	1	Spring guide disc	931118
43	1	Euroseal electronic board 2001 230V	856059
		Euroseal electronic board 2001 110V	856061
44	1	PVC label diam.17 heat symbol	948102
45	1	Painted sheet metal front guard	930188
46	1	Cable 3x0.75 Schuko plug pr.vde 2MT	850040
		Cable AWG 16 American plug pr.vde 2MT	850042
47	1	PVC front panel	948227

SCRAPPING INSTRUCTIONS

The thermosealing machine EUROSEAL® 2001 PLUS has been manufactured using iron materials, electronic components and plastic materials. No particular precautions are required to throw the unit away. Do not dump the unit in an unguarded place; consign it to a waste disposal agency for scrapping. Always refer to the specific laws in the country of use when scrapping equipment.

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